



Cat nr AE00108

Product performance

Product Datasheet

Recombinant Mouse Antibody, clone rIGA/764 to:

CD79a, CD79a molecule

B-cell antigen receptor complex-associated protein alpha chain; CD_antigen: CD79a; Ig-alpha; MB-1 membrane glycoprotein; Membrane-bound immunoglobulin-associated protein; Surface IgM-associated protein; CD79a; IGA; MB1; MB-1


Cellular localization	Plasma membrane, cell surface
Official Symbol (Gene) GeneID SwissProt	CD79A 973 P11912
Confirmed Applications Positive controls Aeonian Rating©	ICC, IHC, PA Daudi, Ramos, Germinal Center B-cells in lymph node, spleen and tonsil. 90
Purification Formulation Amount Isotype Confirmed species reactivity Immunogen Epitope	By Protein G from bioreactor concentrate 200ug IgG/ml in PBS, 0.05% BSA, 0.05% azide (20ug or 100ug) 1mg IgG/ml in PBS (100ug or contact us for quotation) 20ug 100ug Mouse IgG1, kappa Human Recombinant full-length human CD79A protein Unknown
Storage instructions	Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve
Expiration	antibody at 4C for 6-12 months, when kept away from direct sun light. Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.
Warranty	This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any warranty is void. In case of non-conformance, user needs to contact us immediately for replacement or refund.
Liability	This product is for in vitro research use only. Any other applications, such as diagnostics or therapeutics, or in vivo experiments, and the validation of this product therein, are solely at the responsibility of the buyer/user.

see next pages

Product data:

ImmunoHistoChemistry (IHC):


This product was successfully used to stain human spleen sections. Recommended concentration: 0.5-1.5 ug/ml

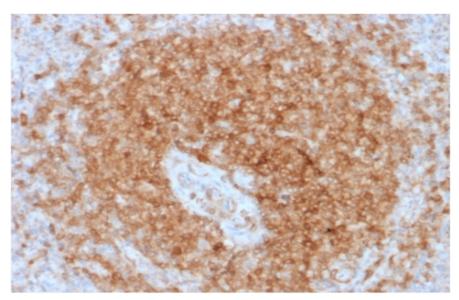

Formaldehyde-fixed, paraffin-embedded human spleen stained with CD79a Recombinant Mouse Antibody AE00108 at 0.5-1.0ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.


ImmunoCytoChemistry (ICC):

This product was successfully used to stain vesicles in Raji. Recommended concentration: 1-3 ug/ml

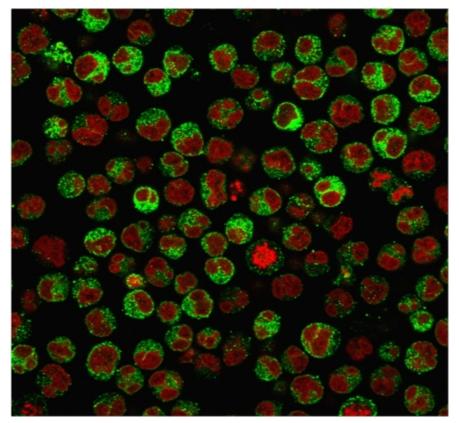


PFA-fixed Raji cells stained with CD79a Mouse Recombinant Antibody AE00108 at 1-2ug/ml for 1h at RT. Detection by confocal microscopy using CF488 (green) for the antibody and RedDot (red) for nuclear staining.

SDS-PAGE Analysis of Purified CD79a Recombinant Mouse Antibody AE00108. Confirmation of Purity and Integrity of Antibody.

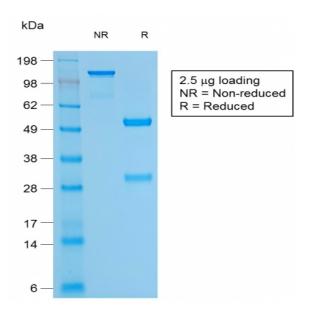

Integrity of the purified antibody AE00108 under non-reduced and reduced conditions, showing intact IgG at around 130kDa (NR) and intact heavy and light chains at 50kDa and 30kDa resp. (R).

Specificity and selectivity of AE00108 to CD79a were tested against >19,000 full-length human proteins on a human protein array. A protein BLAST search against H. sapiens revealed the following closely related other protein: CD79B. This related protein was part of the array used and showed no cross-reactivity signals.

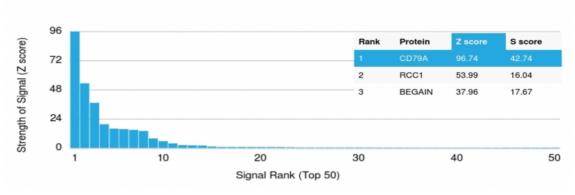

Cross-reactivity assessment of CD79a Mouse Recombinant Antibody AE00108 (1ug/ml) on CDI's Protein Array containing more than 19,000 full-length human proteins.

The Z-score represents the strength of a signal that an antibody (through a fluorophore-tagged secondary reagent) produces when binding to a particular protein on the array. Z-scores are in units of standard deviations (SD's) above the mean value of all signals generated on that array. When Z-scores are arranged in descending order, the difference between two successive values will be the S-score for the first. Thus, the S-score represents the relative specificity of the antibody to its intended target. An antibody is considered specific to its intended target, when it has an S-score of at least 2.5. For example, if an antibody binds to intended protein X with a Z-score of 43 and to the cross-reacting protein Y with a next Z-score of 14, then the S-score for the antibody to intended target X equals 29 (43-14).