



## Cat nr AE00119

Product performance

## **Product Datasheet**

Recombinant Rabbit Antibody, clone MYOD1/2075R to:

## MyoD1, myogenic differentiation 1

Class C basic helix-loop-helix protein 1; myogenic differentiation 1; BHLHC1; MYF3; MYOD; MYOD1; PUM


Cellular localization	Nucleus
Official Symbol (Gene) GeneID SwissProt	MYOD1 4654 P15172
Confirmed Applications Positive controls  Aeonian Rating©	IHC, PA Skeletal muscle, rhabdomyosarcoma 90
Acoman Nating	30
Purification Formulation  Amount Isotype	By Protein A from bioreactor concentrate  200ug IgG/ml in PBS, 0.05% BSA, 0.05% azide (20ug or 100ug)  1mg IgG/ml in PBS (100ug or contact us for quotation)  20ug  Rabbit IgG
Confirmed species reactivity	Human
Immunogen	Recombinant full-length human MyoD1 protein
Epitope	Unknown
Storage instructions	Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C
	or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve antibody at 4C for 6-12 months, when kept away from direct sun light.
Expiration	
Expiration Warranty	antibody at 4C for 6-12 months, when kept away from direct sun light.  Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.  This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions.  User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any warranty is void. In case of non-conformance, user needs to contact us immediately
	antibody at 4C for 6-12 months, when kept away from direct sun light.  Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.  This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any

see next pages

## **Product data:**


ImmunoHistoChemistry (IHC):

This product was successfully used to stain human rhabdomyosarcoma sections. Recommended concentration: 0.3-1ug/ml

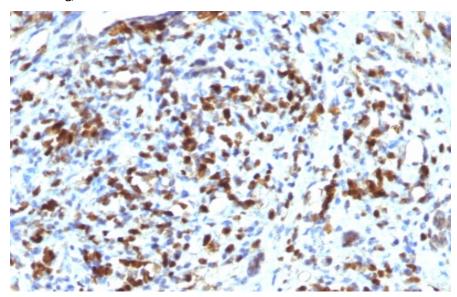


Formaldehyde-fixed, paraffin-embedded human rhabdomyosarcoma stained with MyoD1 Recombinant Rabbit Antibody AE00119 at 0.5-1ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH9 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.

SDS-PAGE Analysis of Purified MyoD1 Recombinant Rabbit Antibody AE00119. Confirmation of Purity and Integrity of Antibody.

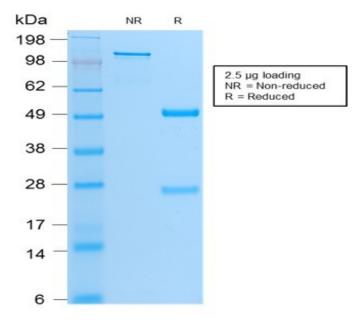

Integrity of the purified antibody AE00119 under non-reduced and reduced conditions, showing intact IgG at around 120kDa (NR) and intact heavy and light chains at 50kDa and 25kDa resp. (R).

Specificity and selectivity of AE00119 to MyoD1 were tested against >19,000 full-length human proteins on a human protein array. A protein BLAST search against H. sapiens revealed the following related other proteins: MYF5 and MYF6. The MYF5 protein was part of the array used and showed no cross-reactivity signals.

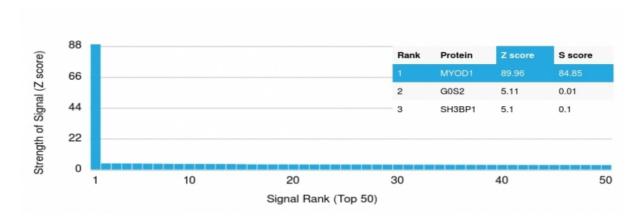

Cross-reactivity assessment of MyoD1 Rabbit Recombinant Antibody AE00119 (1ug/ml) on CDI's Protein Array containing more than 19,000 full-length human proteins.

The Z-score represents the strength of a signal that an antibody (through a fluorophore-tagged secondary reagent) produces when binding to a particular protein on the array. Z-scores are in units of standard deviations (SD's) above the mean value of all signals generated on that array. When Z-scores are arranged in descending order, the difference between two successive values will be the S-score for the first. Thus, the S-score represents the relative specificity of the antibody to its intended target. An antibody is considered specific to its intended target, when it has an S-score of at least 2.5. For example, if an antibody binds to intended protein X with a Z-score of 43 and to the cross-reacting protein Y with a next Z-score of 14, then the S-score for the antibody to intended target X equals 29 (43-14).