



## Cat nr AE00140

Product performance

Recombinant Rabbit Antibody, clone IG507R to:

## pan-IgG, Immunoglobulin G heavy chain

Immunoglobulin heavy constant gamma 1, 2, 3, 4


Cellular localization	Plasma membrane, cell surface and secreted
Official Symbol (Gene) GeneID SwissProt	IGHG 3500, 3501, 3502, 3503 P01857, P01859, P01860, P01861
Confirmed Applications Positive controls  Aeonian Rating©	IHC, PA Raji, PBL, tonsil, spleen
Acoman Nathige	30
Purification Formulation  Amount Isotype	By Protein A from bioreactor concentrate  200ug IgG/ml in PBS, 0.05% BSA, 0.05% azide (20ug or 100ug)  1mg IgG/ml in PBS (100ug or contact us for quotation)  20ug  100ug  Rabbit IgG
Confirmed species reactivity	Human
Immunogen	Recombinant human full-length IGHG protein
Epitope	unknown
Storage instructions	Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve antibody at 4C for 6-12 months, when kept away from direct sun light.
Expiration	Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.
Warranty	This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any warranty is void. In case of non-conformance, user needs to contact us immediately for replacement or refund.
Liability	This product is for in vitro research use only. Any other applications, such as diagnostics or therapeutics, or in vivo experiments, and the validation of this product
	therein, are solely at the responsibility of the buyer/user.

see next pages

## **Product data:**


ImmunoHistoChemistry (IHC):

This product shows select staining of IgG-expressing B cells of human tonsil sections. Recommended concentration: 0.3-1ug/ml

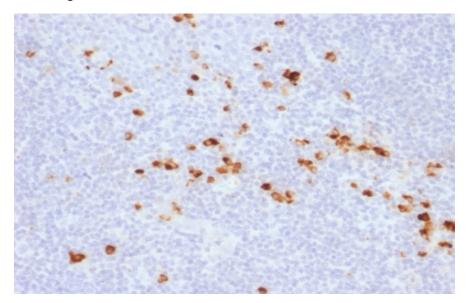


Formaldehyde-fixed, paraffin-embedded human tonsil stained with pan-IgG Rabbit Recombinant Antibody AE00140 at 0.5-1ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.

SDS-PAGE Analysis of Purified pan-IgG Rabbit Recombinant Antibody AE00140. Confirmation of Purity and Integrity of Antibody.

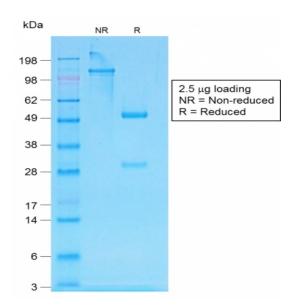

Integrity of the purified antibody AE00140 under non-reduced and reduced conditions, showing intact IgG at around 150kDa (NR) and intact heavy and light chains at 50kDa and 30kDa resp. (R).

Specificity and selectivity of AE00140 to pan-IgG were tested against >19,000 full-length human proteins on a human protein array. A protein BLAST search against H. sapiens reveal no closely related proteins. IGHG2 and IGHG4 were not part of the array, while IGHG3 was.

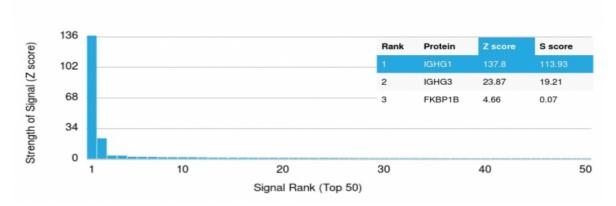

Cross-reactivity assessment of pan-IgG Rabbit Recombinant Antibody AE00140 (1ug/ml) on CDI's Protein Array containing more than 19,000 full-length human proteins.

The Z-score represents the strength of a signal that an antibody (through a fluorophore-tagged secondary reagent) produces when binding to a particular protein on the array. Z-scores are in units of standard deviations (SD's) above the mean value of all signals generated on that array. When Z-scores are arranged in descending order, the difference between two successive values will be the S-score for the first. Thus, the S-score represents the relative specificity of the antibody to its intended target. An antibody is considered specific to its intended target, when it has an S-score of at least 2.5. For example, if an antibody binds to intended protein X with a Z-score of 43 and to the cross-reacting protein Y with a next Z-score of 14, then the S-score for the antibody to intended target X equals 29 (43-14).