



Cat nr AE00169

Mouse Monoclonal Antibody, clone EMD/2168 to:

Emerin

Emerin; EDMD, EMD; LEMD5, STA

Cellular localization Nuclear membranes

Official Symbol (Gene) EMD
 GenelD 2010
 SwissProt P50402

Confirmed Applications IHC, PA, WB
 Positive controls K562, HeLa, Jurkat, kidney, skin.
 Aeonian Rating© 83

Purification By Protein A from bioreactor concentrate
 Formulation 200ug IgG/ml in PBS, 0.05% BSA, 0.05% azide (20ug or 100ug)
 1mg IgG/ml in PBS (100ug or contact us for quotation)
 Amount 20ug 100ug
 Isotype Mouse IgG2b, kappa
 Confirmed species reactivity Human
 Immunogen Recombinant fragment around aa 56-167 of human Emerin protein
 (exact sequence is proprietary)
 Epitope Around aa 56-167

Storage instructions Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve antibody at 4C for 6-12 months, when kept away from direct sun light.

Expiration Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.

Warranty This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any warranty is void. In case of non-conformance, user needs to contact us immediately for replacement or refund.


Liability This product is for in vitro research use only. Any other applications, such as diagnostics or therapeutics, or in vivo experiments, and the validation of this product therein, are solely at the responsibility of the buyer/user.

Product performance see next pages


Product data:

ImmunoHistoChemistry (IHC):


This product shows nuclear membrane staining in sections of breast carcinoma, basal cell carcinoma and renal carcinoma. Recommended concentration: 1-3ug/ml

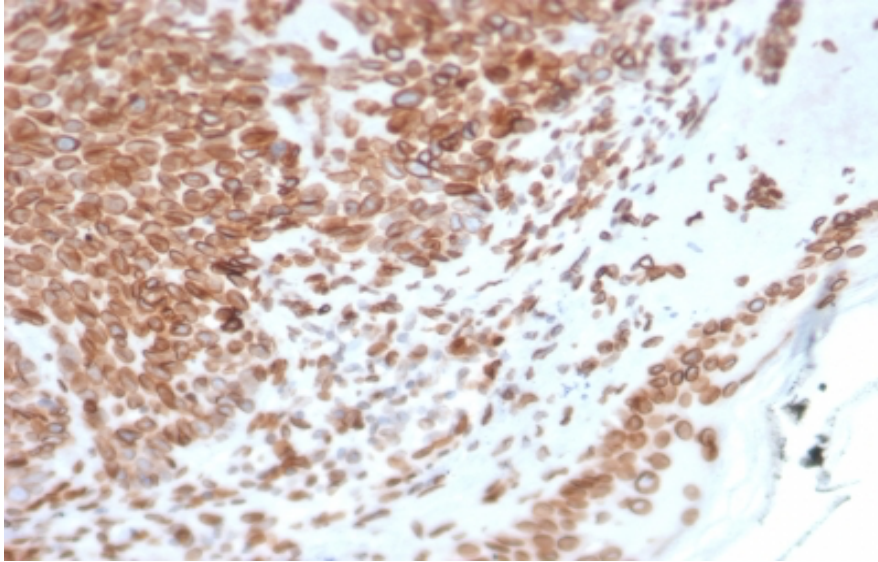

Formaldehyde-fixed, paraffin-embedded human breast carcinoma stained with Emerin Mouse Monoclonal Antibody AE00169 at 1-2ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.

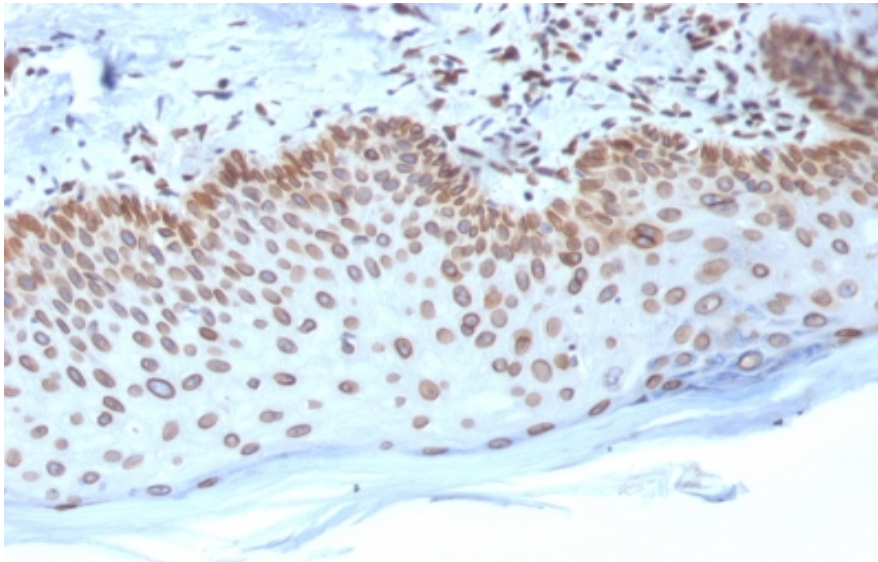


Formaldehyde-fixed, paraffin-embedded human basal cell carcinoma stained with Emerin Mouse Monoclonal Antibody AE00169 at 1-2ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.

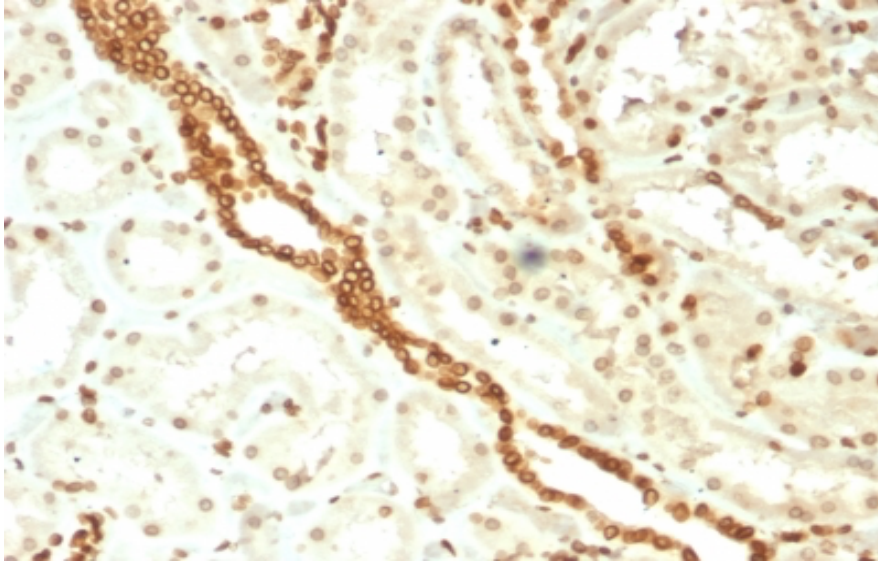

Formaldehyde-fixed, paraffin-embedded human renal carcinoma stained with Emerin Mouse Monoclonal Antibody AE00169 at 1-2ug/ml for 30 minutes at RT. Epitope retrieval: Boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.


Western Blot (WB):

This product was successfully used to stain an approx. 38kDa band in HeLa lysates. Recommended concentration: 0.3-1ug/ml

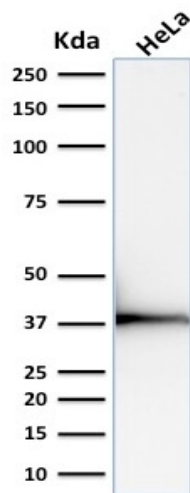


Western Blot of a HeLa lysate (30ug) stained with Emerin Mouse Monoclonal Antibody AE00169 at 0.5ug/ml (1h at ambient temp). ECL staining by HRP.

SDS-PAGE Analysis of Purified Emerin Mouse Monoclonal Antibody AE00169. Confirmation of Purity and Integrity of Antibody.

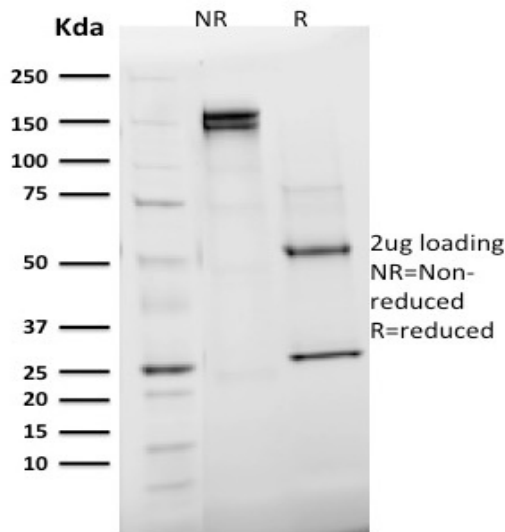

Integrity of the purified antibody AE00169 under non-reduced and reduced conditions, showing intact IgG at around 150+140kDa (NR) and intact heavy and light chains at 50kDa and 25kDa resp. (R). The double band is explained by different states of glycosylation.

Specificity and selectivity of AE00169 to Emerin were tested against >19,000 full-length human proteins on a human protein array. A protein BLAST search against H. sapiens revealed no closely related other proteins.

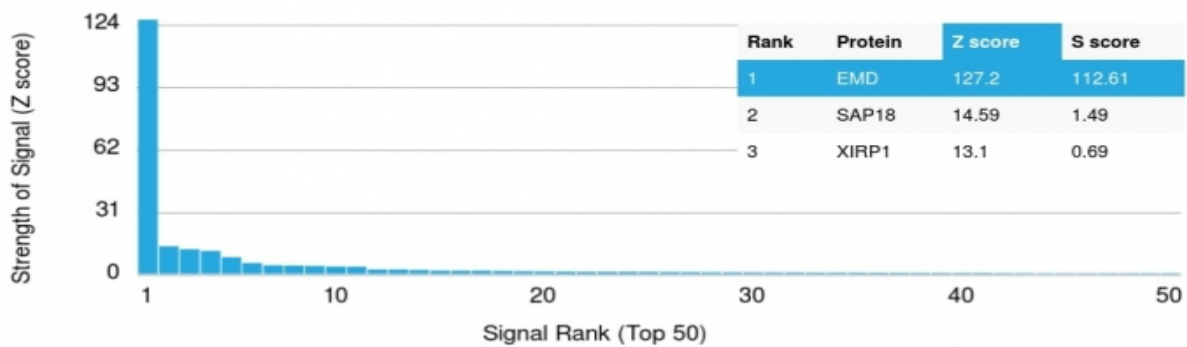

Cross-reactivity assessment of Emerin Mouse Monoclonal Antibody AE00169 (1ug/ml) on CDI's Protein Array containing more than 19,000 full-length human proteins.

The Z-score represents the strength of a signal that an antibody (through a fluorophore-tagged secondary reagent) produces when binding to a particular protein on the array. Z-scores are in units of standard deviations (SD's) above the mean value of all signals generated on that array. When Z-scores are arranged in descending order, the difference between two successive values will be the S-score for the first. Thus, the S-score represents the relative specificity of the antibody to its intended target. An antibody is considered specific to its intended target, when it has an S-score of at least 2.5. For example, if an antibody binds to intended protein X with a Z-score of 43 and to the cross-reacting protein Y with a next Z-score of 14, then the S-score for the antibody to intended target X equals 29 (43-14).