



## Cat nr AE00245

## **Product Datasheet**

Recombinant Rabbit Antibody, clone CHGA/1815R to:

CHGA Chromogranin A; Chromogranin-A; Pituitary secretory protein I; CGA; CHGA; SP-I Precursor to: Vasostatin-1; Vasostatin I; Vasostatin-2; Vasostatin II; EA-92; ES-43; Pancreastatin; SS-18; WA-8; WE-14; LF-19; Catestatin; SL21; AL-11; GV-19; GR-44; ER-37; GE-25; Serpinin-RRG; Serpinin; p-Glu serpinin precursor Cellular localization Extracellular, secreted, secretory vesicles Official Symbol (Gene) CHGA GeneID 1113 SwissProt P10645 **Confirmed Applications** IHC, PA Positive controls pancreas Aeonian Rating© 90 Purification By Protein A from bioreactor concentrate Formulation 200ug lgG/ml in PBS, 0.05% BSA, 0.05% azide (20ug or 100ug) 1mg lgG/ml in PBS (100ug or contact us for quotation) **Amount** 20ug 100ug Rabbit IgG Isotype Confirmed species reactivity Human **Immunogen** Recombinant full-length human CHGA protein **Epitope** Unknown Storage instructions Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve

antibody at 4C for 6-12 months, when kept away from direct sun light.

Expiration Integrity warranted for 24 months after purchase when handled and stored

according to instructions, see below.

Warranty This product is only warranted for the specifications as described in this product

sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any

warranty is void. In case of non-conformance, user needs to contact us immediately

for replacement or refund.

Liability This product is for in vitro research use only. Any other applications, such as

diagnostics or therapeutics, or in vivo experiments, and the validation of this product


therein, are solely at the responsibility of the buyer/user.

Product performance see next pages

## **Product data:**


ImmunoHistoChemistry (IHC):

This product was successfully used to stain Islets of Langerhans in human pancreas sections. Recommended concentration: 1-3ug/ml

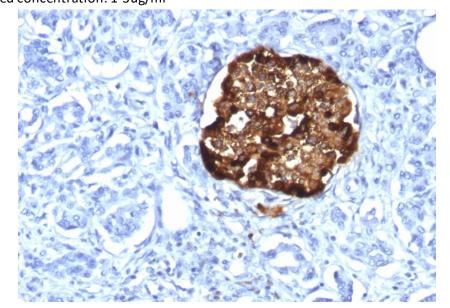


Formaldehyde-fixed, paraffin-embedded human pancreas stained with CHGA Rabbit Recombinant Antibody AE00245 at 1-2ug/ml for 30 minutes at RT. Epitope retrieval: boiling at pH6 for 10-20 min followed by 20 min cooling. DAB staining by HRP polymer.

## SDS-PAGE Analysis of Purified CHGA Recombinant Rabbit Antibody AE00245. Confirmation of Purity and Integrity of Antibody.

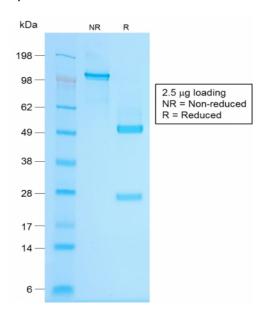

Integrity of the purified antibody AE00245 under non-reduced and reduced conditions, showing intact IgG at around 120kDa (NR) and intact heavy and light chains at 50kDa and 25kDa resp. (R).

Specificity and selectivity of AE00245 to CHGA were tested against >19,000 full-length human proteins on a human protein array. A protein BLAST search against H. sapiens revealed no closely related other proteins.

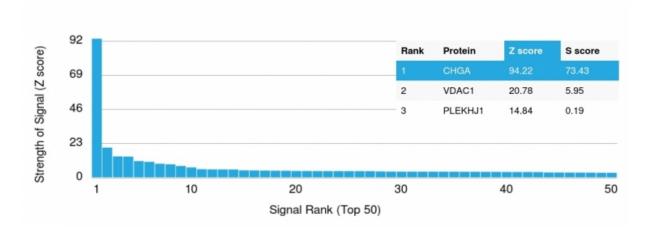

Cross-reactivity assessment of CHGA Rabbit Recombinant Antibody AE00245 (1ug/ml) on CDI's Protein Array containing more than 19,000 full-length human proteins.

The Z-score represents the strength of a signal that an antibody (through a fluorophore-tagged secondary reagent) produces when binding to a particular protein on the array. Z-scores are in units of standard deviations (SD's) above the mean value of all signals generated on that array. When Z-scores are arranged in descending order, the difference between two successive values will be the S-score for the first. Thus, the S-score represents the relative specificity of the antibody to its intended target. An antibody is considered specific to its intended target, when it has an S-score of at least 2.5. For example, if an antibody binds to intended protein X with a Z-score of 43 and to the cross-reacting protein Y with a next Z-score of 14, then the S-score for the antibody to intended target X equals 29 (43-14).