



Cat nr AE00307

Product Datasheet

Rabbit Polyclonal Phospho-Specific Antibody, EGFR-pY1086 to:

EGFR, phosphorylated Y1086


Epithelial Growth Factor Receptor; EGFR; ERBB; ERBB1; HER1; NISBD2; PIG61; mENA

Cellular localization Official Symbol (Gene) GeneID SwissProt	Plasma membrane, cell surface. Translocation to nucleus upon EGF stimulation. EGFR 1956 P00533
Confirmed Applications Positive controls Customers positive feedback Aeonian Rating®	ICC, WB Rat liver, LPS-treated HepG2 70
Purification	Antigen-affinity purified from serum by sequential phosphopeptide and non-phospho-peptide affinity columns.
Formulation	1 mg lgG/ml in PBS, 0.02% azide, 50% glycerol
Amount Isotype Confirmed species reactivity Immunogen	abbit IgG Human, mouse, rat A synthetic peptide around phosphorylated Y1110 of human EGFR (NP_005219.2) corresponding to Y1086 of the mature protein (starts with L25).
Epitope	A short sequence around phosphorylated Y1086 (exact sequence is proprietary)
Storage instructions	Avoid repeated freeze/thaw cycles. For long term storage, keep small aliquots at -20C or -80C and keep one aliquot at 4C for daily experimentations. Azide will preserve antibody at 4C for 6-12 months, when kept away from direct sun light.
Expiration	Integrity warranted for 24 months after purchase when handled and stored according to instructions, see below.
Warranty	This product is only warranted for the specifications as described in this product sheet and only when the product is handled and stored according to instructions. User should validate this antibody in the application and tissue/cell type as required, after confirmation of integrity upon receipt is obtained by reproducing the performance as described below. Should such confirmation not be attempted, any warranty is void. In case of non-conformance, user needs to contact us immediately for replacement or refund.
Liability	This product is for in vitro research use only. Any other applications, such as diagnostics or therapeutics, or in vivo experiments, and the validation of this product therein, are solely at the responsibility of the buyer/user.

Product data:

Western Blot (WB):


This product was successfully used to stain an approx. 160kDa band in mouse and rat liver lysates. Recommended concentration: 1-3ug/ml

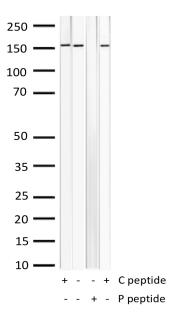

Western Blot of 10ug liver lysates from mouse (lanes 1 and 2) and rat (lanes 3 and 4) stained with EGFR-pY1086 Rabbit Polyclonal Antibody AE00307 at 1ug/ml (overnight at $4^{\circ}C$) in the presence (+) or absence (-) of its antigen (P peptide) or its non-phosphorylated control (C peptide). ECL staining by HRP.

ImmunoCytoChemistry (ICC):

This product shows nuclear membrane staining in LPS-treated cell line HepG2. Recommended concentration: 5-15ug/ml

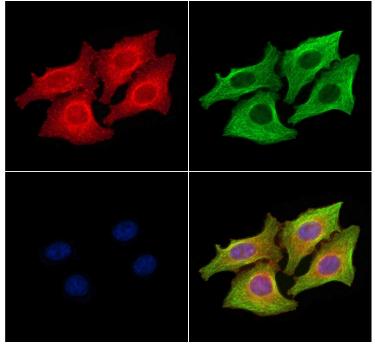

Confocal microscopy of PFA-fixed, and LPS-treated for 4h, HepG2 (0.1% TX100 permeabilised) stained with EGFR-pY1086 Rabbit Polyclonal Antibody AE00307 at 5ug/ml (1h at 37°C). Detection by AF594 (red) for the EGFR-pY1086 rabbit antibody, AF488 (green) for beta-tubulin mouse antibody, and DAPI (blue) for nuclear staining.